



# K9 HERO

Canine Medical Trainer.



**TACTICAL MEDICAL  
SOLUTIONS®**

PHONE 864.224.0081  
TOLL FREE 1.888.822.6331  
FAX 864.224.0064  
E-MAIL [info@tacmedsolutions.com](mailto:info@tacmedsolutions.com)  
WEB [tacmedsolutions.com](http://tacmedsolutions.com)

1250 Harris Bridge Road, Anderson, SC 29621




Tactical Medical Solutions® is a global leading developer and manufacturer of pre-hospital medical gear for the tactical and civilian world.

We currently have over 2,000 products available and a presence in 70 countries around the world. Here at home, we service 7 of the top 10 police departments, including the NYPD and LAPD, as well as major fire and EMS agencies.

We're proud to say that when customers think of TacMed™, they think uncompromising - uncompromising commitment to delivering top-quality products, tactical knowledge, customized training, and reliability in the field.

# K9 HERO

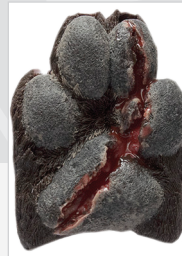

The TraumaFX® Canine Medical Trainer (K9 HERO) is a full-body medical trainer for Operational Canine First Responders, Military Working Dog (WMD) handlers, Veterinarians, and veterinary technicians. Designed in partnership with the Department of Defense (DoD), K9 HERO is a state of the art skills trainer for novices and experienced alike. All training sites are realistic to replicate the look, feel, and function of actual medical procedures. The head of the simulator has been hand-sculpted to mimic a Belgian Malinois, and the exterior is covered in durable, non-staining synthetic fur. K9 HERO allows learners to perform critical life-saving tasks such as maintaining an airway, needle decompression/thoracocentesis, hemostasis, IV insertion, Intraosseous (I/O) infusion, Cardiopulmonary Resuscitation (CPR), tracheostomy, and bandaging.

PRODUCT NUMBER	KGS-TFX-K9-HERO-1	
NSN #	6509016597165	
DIMENSIONS	Length 42" (excludes tail)	
WEIGHT	FULL	50 lbs
	EMPTY	48 lbs
POWER SUPPLY	One (1) 12V Li-ion battery	
BLOOD RESERVOIR	-500 ml	
CASE DIMENSIONS	51"x26"x12"; 115 lbs	

## CONTRACT VEHICLES

GSA	GS-07F-063DA
DLA ECAT	SPE2DH-18-D-0008
PEO STRI TATT II	W900KK-19-D-0005
PEO STRI VPSS	W900KK-18-D-0012
NATO Logistics Stock Exchange	LB-KGS-01N



## K9 HERO KEY BENEFITS

<ul style="list-style-type: none"> <li>▶ Lifelike, full anatomical skeletal motion</li> <li>▶ 8 Intervention training sites</li> <li>▶ IV training site with flash cue</li> <li>▶ Needle decompression/thoracocentesis site allows withdrawal of either air, fluid, or both</li> <li>▶ Adjustable femoral pulses via remote control</li> <li>▶ Reusable lacerated paw for suture and bandage training</li> <li>▶ Realistic airway with biomimetic jaw opening and tongue for proper intubation and bag ventilation training (accepts 10-12mm tubes)</li> <li>▶ Chest rise and fall with correct intubation (stomach rise and fall if incorrectly placed in esophagus)</li> <li>▶ CPR chest compressions with resistance</li> </ul>	<ul style="list-style-type: none"> <li>▶ Optional mouth wound kit with injured mouth and tongue with throat blockage</li> <li>▶ Mouth-to-snout resuscitation</li> <li>▶ Infusible tibial I/O training site with reusable bone plug</li> <li>▶ Saphenous wound requiring appropriate pressure to occlude bleeding (success/failure immediately displayed on the remote control)</li> <li>▶ Tracheostomy training site with multi-layered tissue</li> <li>▶ Realistic look and feel with a weight of ~50 lbs</li> <li>▶ Hydrophobic fur with washable leg socks</li> <li>▶ Microchip insertion site</li> <li>▶ Optional floppy ears</li> </ul>
--	--

# K9 HERO

## K9 HERO FEATURES GUIDE



WIRELESS HANDHELD  
REMOTE (RC) WITH  
PROGRAMMABLE  
SCENARIOS

\* Optional wound kit  
sold separately.



- Accurate oral and nasal airway
- Tracheostomy
- Cephalic IV
- Lacerated front paw (right)
- Needle decompression/  
thoracocentesis site
- Ear tattoo
- Saphenous wound requiring  
appropriate pressure to occlude  
bleeding
- Lifelike, full skeletal motion
- CPR-chest compressions
- Microchip insertion site
- Tibial I/O
- Bilateral femoral pulses

All TraumaFX® products are  
handcrafted in the USA.

**Confidentiality Notice:**

TraumaFX® and Multiple Amputation Trauma Trainer (MATT)® are registered trademarks of TraumaFX® Solutions, Inc. This document contains protected information and its contents constitute Confidential and Proprietary Information. Any unauthorized use, disclosure or distribution is strictly prohibited without prior written consent by an authorized TraumaFX® associate.